



Working together to inspire your patients



Penlet at a glance

	 Penlet	 i-pen	 i-pen[®]	 Axis-D	 Softpen
					
 Manual injection		●	●	●	
 Auto injection	●				●
 Reusable		●	●		●
 Disposable	●			●	
 Hidden needle	●				●
 Multiple injections	●	●	●	●	●
 Variable dosage		●	●	●	●
 Fixed dosage	●				

Penlet System


The only disposable multi-injection autoinjector for 3 ml cartridges. The Haselmeier Penlet is a fully automated, multi-dose, fixed dose injection device designed for use with a standard 3 ml cartridge. The Penlet features Haselmeier's hidden needle design which hides the needle from the patient before and during injection. Pressing the clip on the pen causes the needle to automatically enter the subcutaneous tissue, followed by delivery of the fixed volume of the solution.




The Penlet incorporates the following features:

- | | |
|--|---|
| ✓ Multiple fixed dose injections | ✓ Easy to use |
| ✓ Hidden needle to overcome the fear of injection | ✓ No dose adjustment required |
| ✓ Automatic needle insertion and injection for people with reduced dexterity | ✓ Economical product with a cost per dose less than a prefilled syringe |


Instructions for the use of Penlet

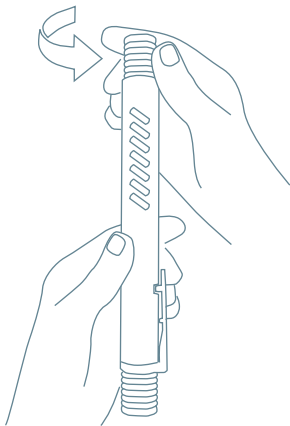


1 / Unscrew pen cap

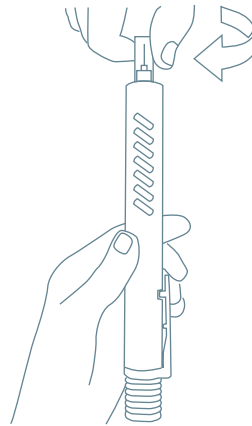


2 / Remove the paper tab and
attach the needle to the pen

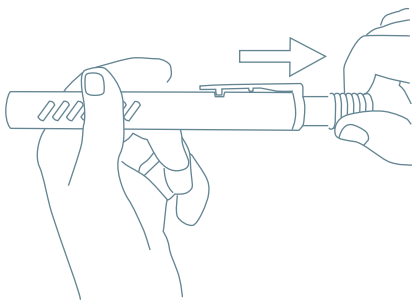


3 / Pull back the plunger until click is heard
and plunger locks into position

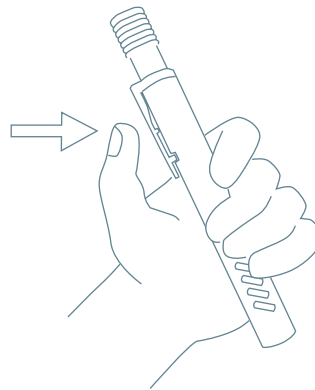

4 / Place the Penlet on injection site
and depress clip

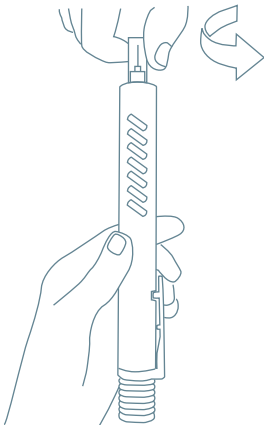

5 / Remove the needle with needle cover

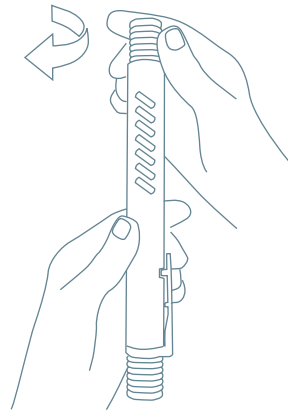

6 / Attach the cap until the next use



Applications

Penlet Application for various
Low-Molecular-Weight Heparins



Pen Features



For injection of Enoxaparin (low-molecular-weight heparins)



Deep vein thrombosis treatment



User Focus

At Haselmeier our mission is to create products enabling a convenient and comfortable experience.

This is why patient feedback is integrated early in our device design.



Bridge to Market

At Haselmeier we provide integrated design, development and industrialization services to help you bridge your serial product into the market.



Manufacturing & LifeCycle Partner

Regardless of your requirements, Haselmeier applies the highest quality standard for manufacturing your drug delivery device to ensure a reliable and reproducible manufacturing and quality process.